


Curapor® transparent Surgical dressing, sterile



Curapor® transparent

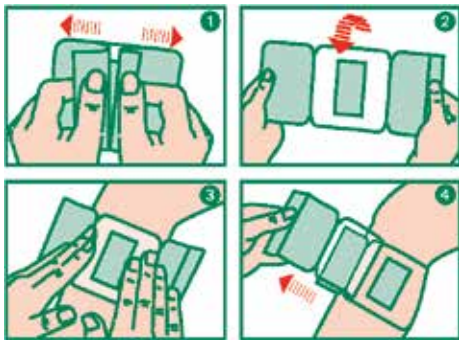

Curapor® Transparent Wound Dressing is easy to apply

- A** just open the backing paper and pull apart ...
- B** ... turn product over ...
- C** ... gently press onto the cleansed skin ...
- D** ... pull off the supportive film – finished!

Fields of application

- sterile treatment of postoperative wounds, injuries sustained in accidents, as well as minor cuts and abrasions
- minimally invasive surgery
- protection against secondary infection
- to cover wounds when showering or swimming

Properties

- ready for use
- transparent and waterproof
- breathable
- impermeable to bacteria
- resistant to common disinfectants
- safe adhesion, precise fit
- leaves little residue; virtually painless removal
- low risk of adhesion to the wound
- also suitable for sensitive skin

Product composition

Carrier film: 100% polyurethane, skin-friendly polyacrylate adhesive (free of colophony and colophony derivatives)
Wound pad: absorbent layer made of viscose/polyester, polyester/polyethylene wound contact layer

To note

Sizes 7 x 5 cm and 10 x 8 cm have rounded edges

Curapor® transparent sterile, individually sealed

Dimensions (cm)	Wound pad (cm)	REF	Shipping units (piece)
7 x 5	4 x 2.5	13099	5/120
10 x 8	6 x 4	13100	5/170
7 x 5	4 x 2.5	13101	50/500
10 x 8	6 x 4	13102	25/200
15 x 10	11 x 5	13103	25/225
20 x 10	16 x 5	13104	25/225
25 x 10	20 x 5	13105	25/225

Not available in all countries